

September 25, 2026



MONTANA TREATMENT COURT CONFERENCE

PharmChek® by SCRAM Systems
Certification for Rural Practitioners





TODAY'S AGENDA:

- ◆ Introductions
- ◆ Completed Training Modules
- ◆ Online Certification
- ◆ General Review:
PharmChek® Sweat Patch
- ◆ Q&A



Products ▾

Resources ▾

Training ▾

About Us ▾

PharmChek Academy

Training Modules

Online Certification



GRAVY

REQUIRED

REQUIRED

Training Modules
Online Certification

Stronger Together:

PharmChek® Officially Joins SCRAM Systems

Learn More 



PharmChek® Products

Learn more about the better, cleaner, smarter option for continuous drug monitoring.

Learn More 



PharmChek® Academy

Comprehensive, virtual, and engaging training sessions to get the most from PharmChek®.

Sign Up 



Online Certification

Take our 25-question exam to certify you have to apply and remove the PharmChek® Sweat Patch

Get Started 

PCA-General Overview
PCA-Application & Removal
PCA-Tampers & Allergies
PCA-Navigating CRL & Reports
PCA-Open Q&A's
PCA-Targeted Roundtable Q&A's
PCA-Guest Speakers

TODAY'S TEST DAY?

I LOVE TESTS!!



Customer Toolbox Pay Your Bill Re-Order PharmChek® Products Request A Quote

PHARMCHEK
by SCRAM Systems

Products ▾ Resources ▾ **Training ▾** About Us ▾

PharmChek Academy
Training Modules
Online Certification

Stronger Together:
**PharmChek® Officially Joins
SCRAM Systems**

Learn More

**PharmChek®
Products**
Learn more about the
better, cleaner, smarter
option for continuous
drug monitoring.
[Learn More](#)

**PharmChek®
Academy**
Comprehensive, virtual,
and engaging training
sessions to get the most
from PharmChek®.
[Sign Up](#)

**Online
Certification**
Take our 25-question
exam to certify you have
to apply and remove the
PharmChek® Sweat Patch
[Get Started](#)



TODAY'S GOAL

**COMPLETE THE REQUIRED
CERTIFICATION FOR
ANYONE THAT IS
APPLYING OR REMOVING
THE
PHARMCHEK® SWEAT
PATCH**

www.pharmchek.com

[Customer Toolbox](#) [Pay Your Bill](#) [Re-Order PharmChek® Products](#) [Request A Quote](#)



[Products](#) [Resources](#) [Training](#) [About Us](#)



PharmChek Sweat Patch Certification Exam

Re-print Certificate of Completion ([click here if you have already passed the exam](#))

To certify that you have been trained to apply and remove the PharmChek® Drugs of Abuse Patch, please fill out all information and take the exam.

You must get at least 22 out of the 25 questions correct to receive a passing score.

First Name *

Last Name *

Agency/Institution *

Customer Type *

Please select best description

Address *

City *

State *

Select a State

Zip Code *

Country *

Select A Country

Work Phone *

Email *

Please only use lowercase letters.

Trained By

PharmChek Staff

Training

PharmChek Academy®

Training Modules +

Online Certification

Customer Support

Do you have a question about your existing account?

Call 817-591-4100
(Then Select Option 2)

You can also contact customer support using the button below.

[Contact Customer Support](#)

PharmChek® Technical Guide

View an online PDF of the PharmChek® Technical Guide

[View Guide](#)

PharmChek® FAQs

View an online PDF of the PharmChek® Sweat Patch Comprehensive FAQs

[View FAQs](#)

REMINDER:

It's an open book test!



TECHNICAL GUIDE

FAQ'S

You must get at least 22 out of the 25 questions correct to receive a passing score.

ONLINE CERTIFICATION



CERTIFICATE *of* PARTICIPATION

John Smith

Has successfully passed the PharmChek® Sweat Patch Drugs of Abuse Training and Certification Course.

Acme Incorporated



Presented By 
PharmChem President & CEO, Kerri Wagner


On This Day March 23, 2022

Advantages of the PharmChek® Sweat Patch: Adults & Juveniles

- ***Extended Detection Window-Continuous Collection***
- **24/7 Accountability**
- **FDA Cleared**
- **Tamper Evident**
- ***No False Positives***
- **Cost-Effective**
- **Removes Testing Gaps**
- ***Confirmed Results (LC-MS/MS)***
- **Parent Drug & Metabolite**
- **Visual Deterrent**
- **Sanction or Reward**
- **Court Validated**
- **35+ Years of Peer-Reviewed Research**
- **Easy Application & Removal**
- **Non-Invasive & Dignified**
- **Hygienic**
- **7–10+ Day Wear**
- **Forensically Reliable & Defensible Results**
- ***Informed Intervention Strategies***

PHARMCHEK® SWEAT PATCH

CONSTRUCTION

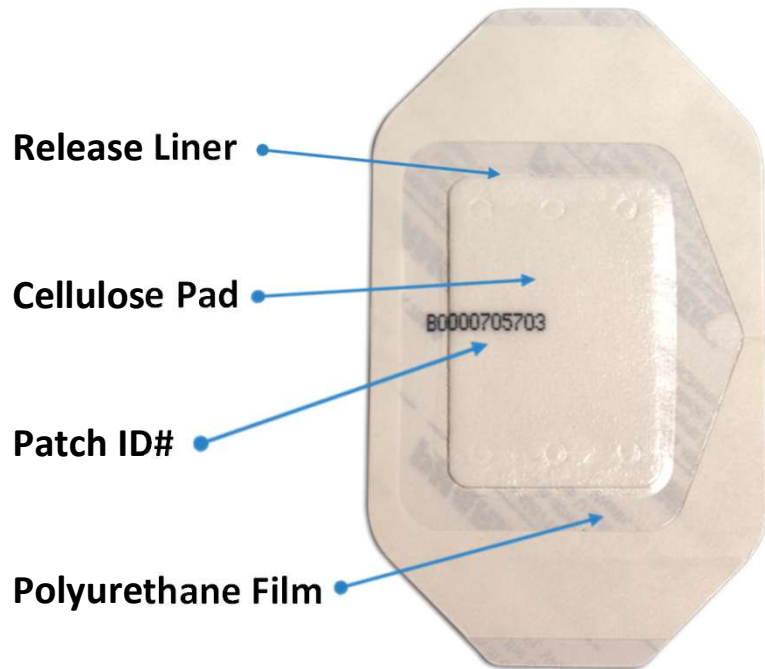

- Transdermal & Continuous
- Collecting insensible perspiration
- Drugs are trapped in the pad
- Semi-permeable membrane allows water vapor & gases to escape
- Skin below remains healthy & sterile
- No evidence of environmental contamination

WHERE & WEAR TIME

The PharmChek® Sweat Patch is worn on the upper arm (lower midriff and lower back are alternative sites).


The Sweat Patch should be worn for a minimum of 24 hours to ensure an adequate amount of sweat is collected for testing.

Recommended wear time is 7-10 days with **SOME** individuals being able to wear it up to 14 days.




PHARMCHEK[®]

EXPANDED PANEL

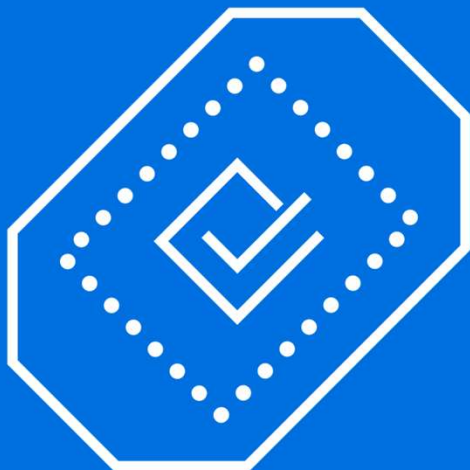

Expanded Drug Panel	LC-MS/MS Confirmation Cutoffs
Methamphetamine AMP/Metabolite	10 ng/mL Positive 2.0-9.9 ng/mL Present
Amphetamine	10 ng/mL
Cocaine BE/Metabolite	10 ng/mL Positive 2.0-9.9 ng/mL Present
Benzodiazepines (9 analogs)	1 ng/mL
Buprenorphine, Norbuprenorphine	0.5 ng/mL
Fentanyl, Nor fentanyl	1 ng/mL
Methadone	5 ng/mL
Codeine	10 ng/mL
Morphine	10 ng/mL
6-AM	10 ng/mL
Hydrocodone, Hydromorphone	10 ng/mL
Oxycodone, Oxymorphone	10 ng/mL
PCP	7.5 ng/mL
THC (Delta-9)	0.5 ng/mL

PHARMCHEK[®] STANDARD PANEL

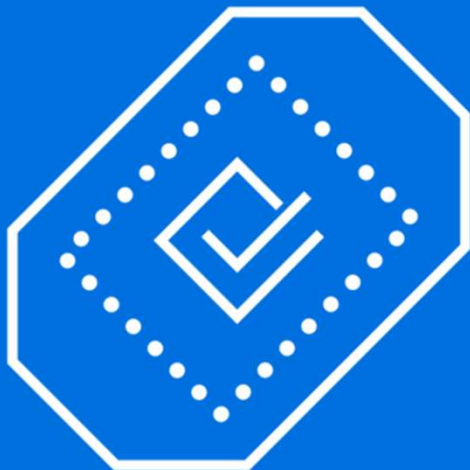

Standard Drug Panel	LC-MS/MS Confirmation Cutoffs
Methamphetamine AMP/Metabolite	10 ng/mL Positive 2.0-9.9 ng/mL Present
Amphetamine	10 ng/mL
Cocaine BE/Metabolite	10 ng/mL Positive 2.0-9.9 ng/mL Present
Codeine	10 ng/mL
Morphine	10 ng/mL
6-AM	10 ng/mL
Hydrocodone, Hydromorphone	10 ng/mL
PCP	7.5 ng/mL
THC (Delta-9)	0.5 ng/mL
Fentanyl Add-on	LC-MS/MS Confirmation Cutoffs
Fentanyl, Nor Fentanyl	1 ng/mL



WHAT IS A CUT-OFF LEVEL?



One of the key concepts within drug testing is the application of a cut-off level.

This is the point which segregates a test result as being either positive or negative. For drug screening tests, a cut-off is chosen that will optimize drug detection but minimize the number of false positive results.

****It is extremely important to note that a negative result in drug testing does not mean NO drugs are present. Rather, NO drugs are present at the “cut-off level”***




CHOOSING THE RIGHT METHOD



Category	PharmChek® Patch	Urine Testing	Oral Fluid Testing	Hair Testing
Key Strengths	• Dual analyte required • 7-10+ day wear • FDA-cleared • LC-MS/MS confirmed • 35 yrs. of scientific research & court acceptance • Dedicated training & customer support • Only pay for labs <u>actually received</u>	• Court-accepted and widely understood • Rapid test cups available • Detect recent use effectively (1-3 days) • Broad drug panel options • Good for spot-in-time testing	• Fast, simple, non-invasive • Rapid test kits available • Gender-neutral • Detects very recent use • Good for spot-in-time testing	• Robust history and peer reviewed studies • Long detection window (mos./yrs) • Ideal for detecting chronic drug use
Key Limitations	• No instant results • Lab shipping not included • Higher cost than rapid UA/OF tests • Smaller panel selection than lab confirmed UA or OF	• Risk of tampering • Not trauma-informed • Shy bladder issues • Short detection window	• Very short detection window • Requires second device for confirmations • Many rapid tests fail to meet performance standards (IFU)	• Blind to recent drug use (10-15 days) • Influenced by hair treatments • Prone to environmental contamination • Requires lab analysis

GAPS IN DETECTION

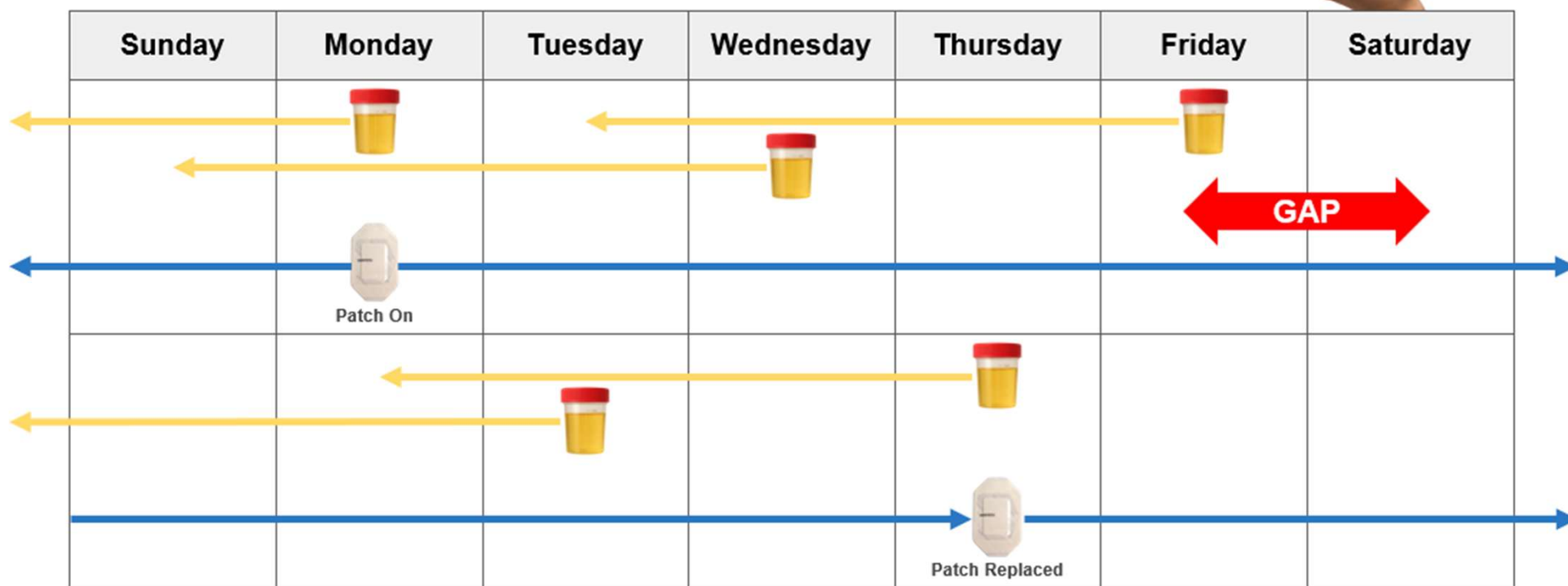

PharmChek® is the only method that tests “forward”

PHARMCHEK[®] LAB TESTING WORKFLOW



ACCESSIONING

Patch sample logged, inspected, barcoded and chain of custody verified



IMMUNOASSAY SCREEN

Drugs extracted and reviewed, presumptive positives move to certification and confirmation, negatives are released by Certifying Scientist



CONFIRMATION

LC-MS/MS confirmation identifies & quantifies drugs and metabolites



CERTIFICATION

Certifying Scientist reviews screening & confirmation results, conducts QA checks, discrepancies are retested or cancelled

CRL DOCUMENTS

CHAIN of CUSTODY

05-Nov-2025 Clinical Reference Laboratory 15:50
CLIA #17D2005163, SAMHSA #0007, CAP #30211-03

PHARMCHEM INC NAME: SAMPLE ID: PB37
DOB: N/S COLLECTED: 10/20/25
ID: 0824 RECEIVED: 11/05/25
GENDER: N/S REPORTED: 11/05/25
SLIP ID: X0003 FAX: (8
REF ID: 75 PANEL ID: WC82
REFERENCE 1:
REFERENCE 2:

PH: (8
COLL. SITE ID: N/S

REASON FOR TESTING: OTHER
SAMPLE TYPE: SWEAT TESTING

INITIAL TEST	RESULT / STATUS	CUTOFF/EXPECTED VALUES
AMPHETAMINE, SWEAT	 NEGATIVE	10 ng/mL
METHAMPHETAMINE, SWEAT	 NEGATIVE	10 ng/mL
BENZODIAZEPINES, SWEAT	 NEGATIVE	2.0 ng/mL
BUPRENORPHINE, SWEAT	 NEGATIVE	1.0 ng/mL
COCAINE METABOLITE, SWEAT	 NEGATIVE	10 ng/mL
PCP, SWEAT	 NEGATIVE	7.5 ng/mL
FENTANYL, SWEAT	 NEGATIVE	1.0 ng/mL
METHADONE, SWEAT	 NEGATIVE	5.0 ng/mL
OPIATES, SWEAT	 NEGATIVE	10 ng/mL
OXYC/OXYM, SWEAT	 NEGATIVE	10 ng/mL
THC, SWEAT	 NEGATIVE	0.8 ng/mL

LAB DIRECTOR: BRITTANY SCOTT, PHD, D-ABFT-FD

REPORT CERTIFIED BY A12

LAB REPORT

PHARMCHEK
Chain of Custody for Analysis of
PharmChek® Drugs of Abuse Sweat Patch
591-4100 | pharmchek.com | 2411 E Loop 820 N | Fort Worth, TX 76118

1 Results Name and Address
Premier Drug Testing
123 Main Street
Small Town, USA 12345

2 PharmChek® Application
FIRST NAME: JOHN
LAST NAME: SMITH
Donor ID: 42050
Observer Name: SARAH JOHNSON
PharmChek® No: 10-20-25
Observer's Initials: SJ
Donor's Initials: JS

3 PharmChek® No
10-20-25

4 Tests Ordered
☒ Standard Panel (WA07) OR ☐ Expanded Panel (WC82)
☐ Fentanyl Add-on (S229)

5 Treatment Status / Reason for Specimen
☒ 01 Random ☐ 02 Probation ☐ 03 Retest
☐ 04 Medical ☐ 05 In Treatment ☐ 06 Pre-Trial
☐ 07 Surveillance (No Treatment)
☐ 08 Other:

6 Observer's Certification at PharmChek® Application
I certify that I applied the PharmChek® identified by the PharmChek® number on this form in accordance with the required procedures.
Observer's Signature: Sarah Johnson

7 Specimen Bag
17
18
19
20
21
SECURITY SEAL
PharmChem, Inc
Barcode
762179706
762179706
762179706
PHARMCHEK®
Donor's Initials
Observer's Initials Date Collected

8 Account Name
Premier Drug Testing
762179706
Account No. 123456789

9 PharmChek® Removal
Date 10-31-25
Observer's Initials SJ
Donor's Initials JS

10 PharmChek® Use Information
Did the PharmChek® appear to be tampered with or compromised?
☐ YES ☒ NO If YES, how?
Comments: R-ARM → G-V

11 Last Use Questionnaire
Name of medications used during PharmChek® wear:
see case file
Dates Used

12 Donor Completes
Donor Certificate and Consent: I certify that the specimen accompanying this form is my own. Further, I certify that the specimen was sealed with a tamper-proof seal in my presence. The information provided on this form, and on the label, is correct. Also, I consent to the laboratory analysis of the specimen accompanying this form and to the release of the laboratory result, as well as the information recorded on this form, to the organization and/or individual listed on this form.
Donor's Signature: John Smith

13 Observer's Certification at PharmChek® Removal
I certify that I removed the PharmChek® identified by the PharmChek® number on this form in accordance with the required procedures. I certify that I applied the numbered security seal and barcode label to the specimen bag in the Donor's presence. I have verified that the specimen number on the form, the barcode, and security seal are identical.
Observer's Signature: Sarah Johnson

22 Instructions to Observer
RECORD application information on this page and print before applying the patch.
APPLY PharmChek® according to Administrative Procedures.
INITIAL form upon PharmChek® application and sign Observer's Certification.
REMOVE PharmChek® according to Administrative Procedures.
AFFIX security seal and barcode to PharmChek® specimen bag as illustrated above.
INITIAL and date security seal. Mail via US Mail to:
CRL Tox Set UP, PO Box 218991 Kansas City, MO 64121
If shipping via FedEx, the ship to address is CRL, 11711 W. 83rd Terrace, Lenexa, KS 66150

05-Nov-2025

Clinical Reference Laboratory
CLIA #17D2005163, SAMHSA #0007, CAP #30211-03

15:50

PHARMCHEM INC

NAME:
DOB: N/S
ID: 0824
GENDER: N/S
SLIP ID: X0003
REF ID: 75
REFERENCE 1:
REFERENCE 2:

SAMPLE ID: PR37

COLLECTED: 10/20/25
RECEIVED: 11/05/25
REPORTED: 11/05/25
FAX: (8
PANEL ID: WC82

PH: (8
COLL. SITE ID: N/S

REASON FOR TESTING: OTHER
SAMPLE TYPE: SWEAT TESTING

INITIAL TEST	RESULT / STATUS	CUTOFF/EXPECTED VALUES
AMPHETAMINE, SWEAT	 NEGATIVE	10 ng/mL
METHAMPHETAMINE, SWEAT	 NEGATIVE	10 ng/mL
BENZODIAZEPINES, SWEAT	 NEGATIVE	2.0 ng/mL
BUPRENORPHINE, SWEAT	 NEGATIVE	1.0 ng/mL
COCAINE METABOLITE, SWEAT	 NEGATIVE	10 ng/mL
PCP, SWEAT	 NEGATIVE	7.5 ng/mL
FENTANYL, SWEAT	 NEGATIVE	1.0 ng/mL
METHADONE, SWEAT	 NEGATIVE	5.0 ng/mL
OPIATES, SWEAT	 NEGATIVE	10 ng/mL
OXYC/OXYM, SWEAT	 NEGATIVE	10 ng/mL
THC, SWEAT	 NEGATIVE	0.8 ng/mL

LAB DIRECTOR: BRITTANY SCOTT, PHD, D-ABFT-FD

COLLECTED:

**Date of
Application**

RECEIVED:

**Date the Lab
Received**

REPORTED:

**Results reported
out & available**

REMOVAL DATE:

Located on COC

PHARMCHEM INC

NAME:
DOB: N/S
ID: CE431
GENDER: N/S
SLIP ID: X0 03 4341
REF ID: 7 6 00002
REFERENCE 1:
REFERENCE 2:

SAMPLE ID: PB 8347
COLLECTED: 08/22/25
RECEIVED: 09/04/25
REPORTED: 09/06/25
FAX: N/S
PANEL ID: WC82

REASON FOR TESTING: OTHER
SAMPLE TYPE: SWEAT TESTING

CCF CONTAINS COLLECTOR COMMENTS. PLEASE REVIEW NON-FEDERAL CCF IN CRL'S WEBOASIS OR OBTAIN A COPY BY EMAILING CLSCS@CRLCORP.COM.

INITIAL TEST	RESULT / STATUS	CUTOFF/EXPECTED VALUES
AMPHETAMINE, SWEAT	NEGATIVE	10 ng/mL
METHAMPHETAMINE, SWEAT	POSITIVE	10 ng/mL
BENZODIAZEPINES, SWEAT	NEGATIVE	2.0 ng/mL
BUPRENORPHINE, SWEAT	NEGATIVE	1.0 ng/mL
COCAINE METABOLITE, SWEAT	NEGATIVE	10 ng/mL
PCP, SWEAT	NEGATIVE	7.5 ng/mL
FENTANYL, SWEAT	NEGATIVE	1.0 ng/mL
METHADONE, SWEAT	NEGATIVE	5.0 ng/mL
OPIATES, SWEAT	NEGATIVE	10 ng/mL
OXYC/OXYM, SWEAT	NEGATIVE	10 ng/mL
THC, SWEAT	NEGATIVE	0.8 ng/mL

CONFIRMATION	RESULT / STATUS	CUTOFF VALUE
MS AMPHETAMINE, SWEAT	PRESENT	
--AMPHETAMINE CUTOFF VALUE--		
2.0 - 9.9 NG/ML	PRESENT	
10.0 NG/ML	POSITIVE	

TEST CERTIFIED BY A90


MS METHAMPH, SWEAT	131.0	POSITIVE	10.0 ng/mL
--------------------	-------	----------	------------

TEST CERTIFIED BY A90

POSITIVE LAB REPORT

INITIAL TEST

CONFIRMATION

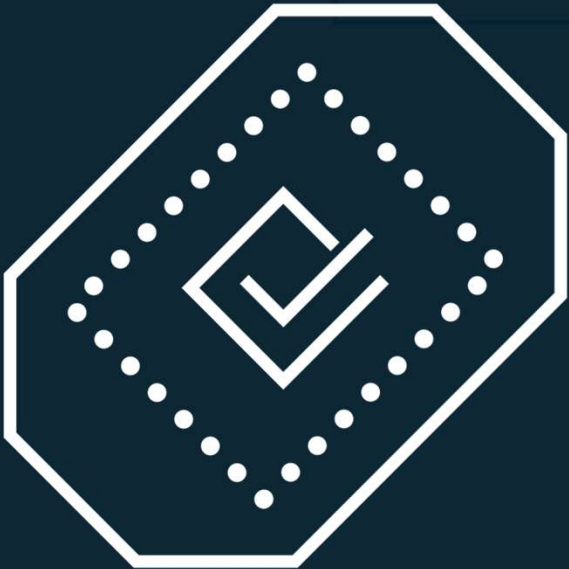

Building client trust, increasing accountability, and improving quality of life are drivers of a successful personalized drug screening program.



DEPTH

(IMPLEMENTING A SWEAT PATCH PROGRAM)

"We selected the successful kids... Out of our first ten tests ... we had five positives for THC, and three for methamphetamine. It was an eye-opener! These weren't individuals we suspected of current use, yet the patch revealed a very different story. PharmChek® was capturing real-time drug use we would have missed with urine testing alone."

Ed Gilligan

Director of Juvenile Probation, Yuma, AZ

D: Do the Math

E: Educate

P: Plan the Roll-out

T: Train

H: Hold the Line

COURT INSTRUCTION & CLIENT RESPONSIBILITY



- Must comply with program rules
- Can still perform normal activities, including work, swimming, showering, exercising, etc.
- Caution: sauna, hot tub, sweat lodge, tanning bed
- If client removes patch without agency/court approved, it's treated as non-compliance
- Nothing can be added to or subtracted from the patch
- Take a picture of the patch when it is applied to compare with when removed



INSTRUCTIONS FOR:

Application

PHARMCHEK
DRUGS OF ABUSE SWEAT PATCH

Always wear gloves when applying or removing the patch.



1 Fill out the application section of the chain-of-custody form (in the donor's presence.) It is mandatory for the observer to sign the chain-of-custody form.



2 Thoroughly cleanse the skin where the patch will be placed with an alcohol pad. The arm should redden as you aggressively cleanse the skin. Continue cleansing with additional pads until the pad is no longer discolored. *If using an overlay, an area at least 5x6 inches must be cleansed with another alcohol wipe to ensure the overlay will adhere properly. Wait 60-90 seconds for the alcohol to dissipate.*



3 Remove the patch from the sealed wrapper. Ensure the arm is flexed while you apply the patch to the donor. Rub your finger around the clear film several times, paying close attention to the corners and edges of the patch.

OPTIONAL OVERLAY



4 Find the slit in the paper border on the outside of the patch. As you peel off the border, press the film to the skin to ensure adhesion.



5 Continue to press the patch to the skin, making sure that the edges and corners are well-adhered to the skin.



6 If using an overlay, remove the full back from the film and place over the top of the patch. Be careful that the only portion showing in the window is the pad. Before removing the quadrant cast liners, rub your index finger around the film for several seconds, paying close attention to ensure the edges and corners are well-adhered.

INSTRUCTIONS FOR:

Removal

PHARMCHEK
DRUGS OF ABUSE SWEAT PATCH

Always wear gloves when applying or removing the patch.



1 Verify the ID# on the outside of the patch matches that of the number written on the donor's original chain-of-custody form.



2 Place your index finger over the corner of the pad, while you instruct the donor to pull down ONLY the top edge of the film to release it from the skin.



3 With the tweezers in hand, continue to pull down the film to expose the top of the pad. Clamp the tweezers over the pad. In one fluid motion, pull up on the pad as you pull down the film.



4 Promptly place ONLY the pad into the specimen bag (the small, clear bag) and seal completely. Discard the film! DO NOT place anything other than the pad into the specimen bag.



Digital CoC Form:

If using a traditional CoC form with a second ply, proceed to step 5B.

Remove a single sticker group from the provided digital chain-of-custody barcode stickers sheet (provided) and apply the sticker group to the space provided on the digital form. Remove one barcode sticker from the chain of custody form and affix it to the space provided at the form's top-right corner.

From this point, follow the remaining instructions beginning in step 5B.




Traditional CoC Form:

Whether you use a digital or traditional CoC form, remove one barcode sticker from the chain of custody form and affix it to the outside of the specimen bag. Complete the chain of custody form in the donor's presence. Both the observer and the donor MUST sign in the required spaces provided.



6 Use the security seal from the chain of custody form to seal the specimen bag. Both the observer and the donor must initial the seal.

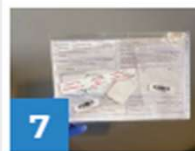

7 Place the specimen bag (containing only the pad) and the second ply of the chain of custody form (if using the traditional CoC form) into the transport bag (the larger of the two bags) and seal completely. DO NOT staple the CoC form to the specimen bag. Affix the final barcode label to the outside and insert it into the brown laboratory mailer (provided).



8 Mail the specimen to one of the three available locations. For timely delivery, the return address must be placed on the top left-hand corner of the shipping envelope.

FedEx GROUND
Clinical Reference Lab, Inc.
11711 W. 83rd Terrace
Lenexa, KS 66214

FedEx EXPRESS
Clinical Reference Lab, Inc.
11711 W. 83rd Terrace
Lenexa, KS 66214

US MAIL
Clinical Reference Lab, Inc.
PO Box 218991
Kansas City, MO 64121

ALL RISE BEST PRACTICE STANDARDS



“**Sweat** testing and other extended detection methods are linked to **higher completion rates** and **lower recidivism**.”

All Rise (2025), Adult Treatment Court Best Practice Standards, p. 29, 87.

“**Sweat** and transdermal alcohol detection devices have demonstrated effectiveness for both **detection** and **deterrence**.”

All Rise (2025), Adult Treatment Court Best Practice Standards, p. 99.



COURT VALIDITY



SAMHSA: Alternative Specimen

“Also, scientific advances in the use of head hair, sweat, and oral fluid in detecting drugs have made it possible for these specimens to be used in Federal programs with the *same level of confidence that has been applied to the use of urine.*”

—SAMHSA, Federal Register, 4/13/04, 69 FR 19689



SAMHSA
Substance Abuse and Mental Health
Services Administration



FDA CLEARANCE



- Cleared by the US FDA under the 510(k) premarket notification process
- Classified as Class II medical device
- Specimen collection for drug testing

IS THE PHARMCHEK® SWEAT PATCH SUBJECT TO PEER-REVIEWED ACADEMIC STUDIES?



Courts regularly cite dozens of peer-reviewed academic studies that support the reliability, admissibility, and forensic defensibility of the PharmChek® Sweat Patch.

Case Study

More than 100 peer-reviewed clinical and scientific papers from research centers worldwide have validated the PharmChek® Sweat Patch as an accurate and reliable drug testing device.



FEDERAL & STATE COURTS



PharmChek® Sweat Patch Court Cases

Search by **STATE** to view Federal, State, or Local court cases involving the PharmChek® Drugs of Abuse Sweat Patch

Filter by State

All States

- ✓ Upheld as a valid and reliable drug testing device in the courts 75+ times
- ✓ PharmChek® Sweat Patch Upheld, U.S. Court of Appeals, 8th Circuit, 2007

For additional information on all court cases, visit pharmchek.com


LITIGATION SUPPORT



Laboratory


- Lab's scope is narrow
- Testify to the result itself
- Testify to lab process
- They DO NOT testify regarding the validity of the collection device

PHARMCHEK
DOES NOT
PROVIDE DIRECT
TESTIMONY

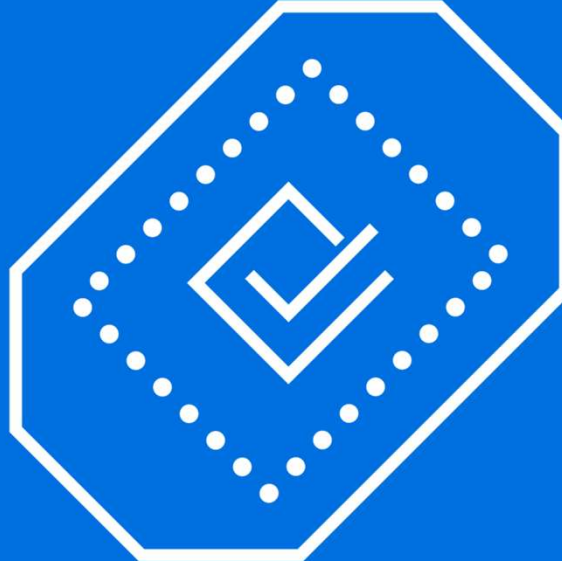

Expert Witness

- Dr. Leo Kadehjian or Dr. Suman Rana
- Reliability of the collection device
- Claims of environmental contamination
- 3rd party: Not associated with CRL or PharmChek

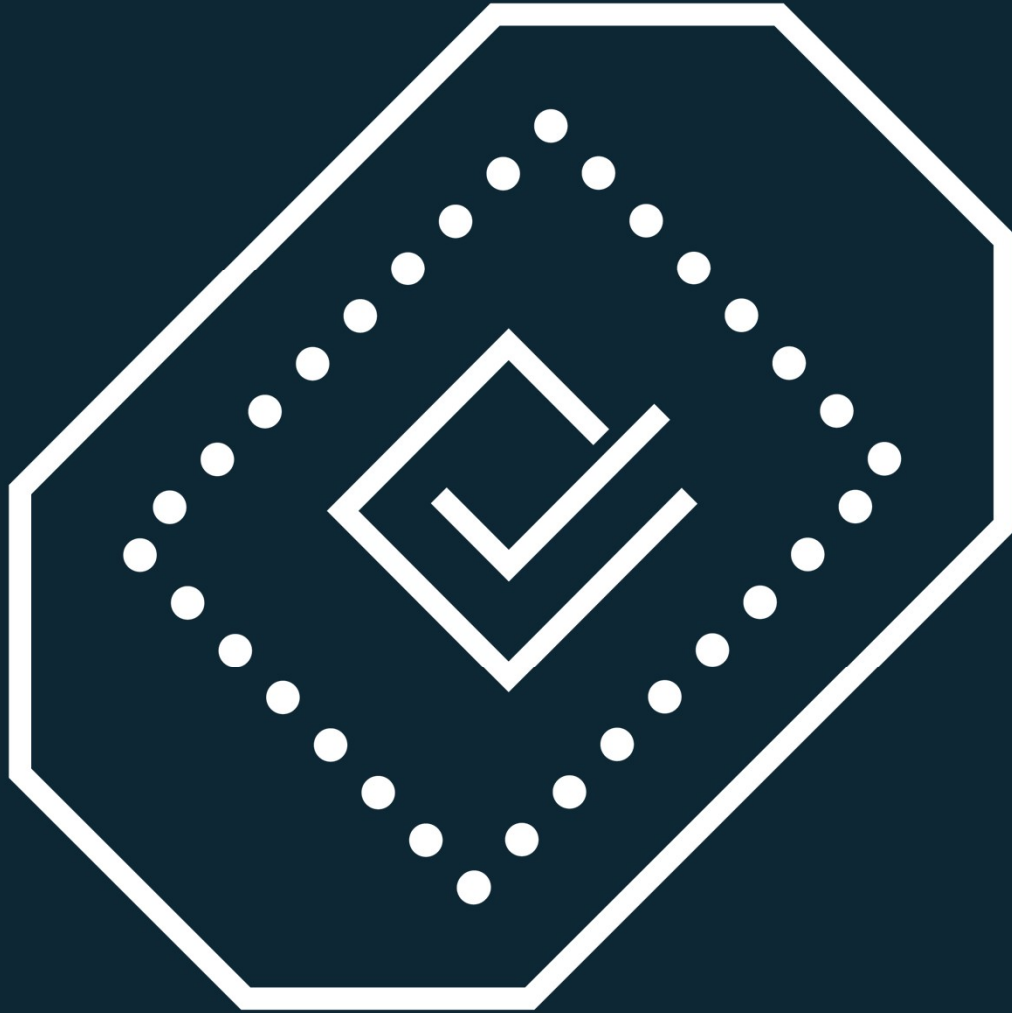

Q & A



Melissa Bosch

SCRAM Account Manager

 **(720) 607-1031**

 **mbosch@scramsystems.com**

Sara McDaniel

SCRAM NW Territory Regional Sales Manager

 **(509) 863-3476**

 **sara.mcdaniel@scramsystems.com**

Kimberly Henderson




PharmChek Director of Customer Success

 **(817) 591-4076**

 **khenderson@scramsystems.com**

pharmchekcustomersupport@scramsystems.com

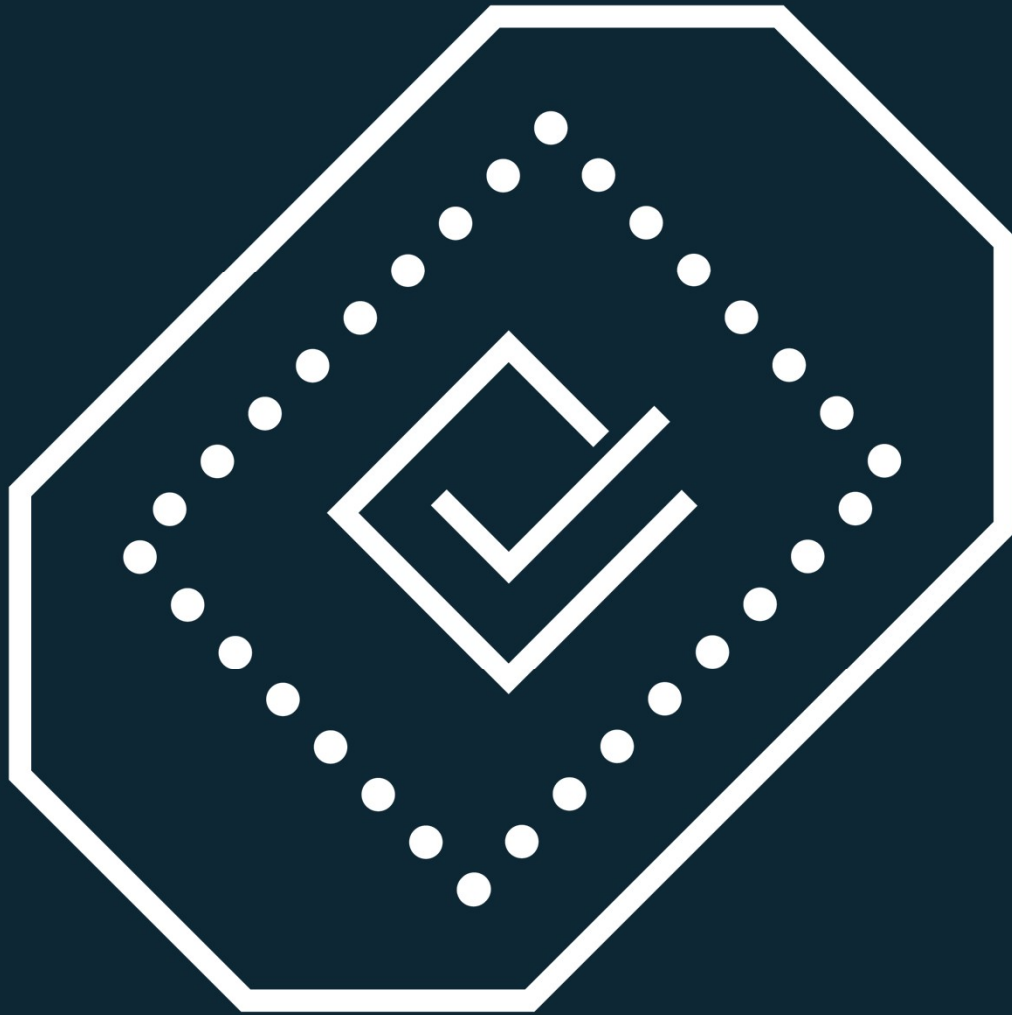

THANK YOU!